


3008 / 3008A Hydraulic Operation Table (Head operating)

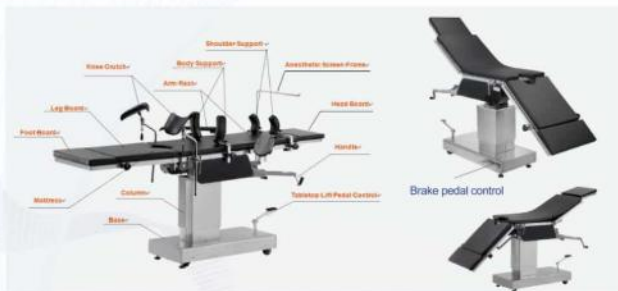

Applications:

The operating table is supplied for thoracic surgery, abdominal surgery, ophthalmology, otorhinolaryngology, obstetrics and gynecology, urinary surgery and orthopaedics.

Features

- Model 3008 is standard type;
Model 3008A: Based on 3008, add base shell, stand column and bed surface made of stainless steel.
 - Head operating universal table.
 - The rise and drop of operating table was handled by pedaling the oil pump, other various movements were handled by the mechanical transmission head.
 - Shoulder support, body support, arm rest, with leather wrap.
- Main Technical Parameter:

Main Technical Parameter:

1. Table Height Range: 760mm~1010mm
2. Length and width: 2010mm x 480mm
3. Lateral Tilt (left/right): $\geq 20^\circ$ / $\geq 20^\circ$
4. Trengelenburg / Reverse: $\geq 25^\circ$ / $\geq 25^\circ$
5. Back Section up/down: $\geq 80^\circ$ / $\geq 10^\circ$
6. Head Section up/down: $\geq 60^\circ$ / $\geq 90^\circ$
7. Leg Section down: $\geq 90^\circ$
9. Waist board lifting: 120 mm

Standard Accessories:

A pair of shoulder support frame, a pair of body support frame, a pair of arm support frame, a pair of leg support frame, a screen frame